

We protect  
your health.



Dr. Schumacher

# Our new player for TEAM ENZYME




## PERFEKTAN® ENZYME

The *PERFEKT* addition for your safety



bactericidal,  
yeastocidal,  
enveloped viruses  
in 5 min.

PREVENT

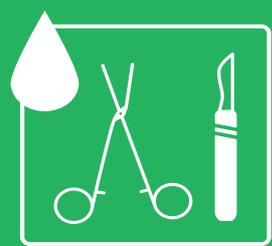

### High level of material compatibility through a neutral pH

The manual reprocessing of sensitive medical products places increased demands on your safety - we have designed PERFEKTAN ENZYME for these situations: In addition to good cleaning performance, you can rely on the disinfection effect of PERFEKTAN ENZYME for your protection.

### Disinfectant effect with pleasant scent

A QAC was deliberately chosen to be the basis of the disinfectant properties and aldehydes as well as amines were omitted from the formulation. Thus, the bactericidal, yeastocidal effectiveness plus activity against enveloped viruses of PERFEKTAN ENZYME is achieved without protein fixation and with very good compatibility to other process chemicals. PERFEKTAN ENZYME can be used in an immersion bath procedure as well as in ultrasonic baths, irrespective of this the discreet scent makes the working conditions more pleasant.

### For effective non-fixing enzymatic pre-cleaning

The strong enzyme complex is not only composed of amylase, lipase and protease, but also contains an additional cellulase. Furthermore, surfactants in this pH-neutral formulation ensure thorough and gentle cleaning with good material compatibility for classical surgical instruments as well as sensitive medical products.

Single Unit	Unit	Content	REF
dispenser bottle	12	1 L	00-128-010E
bottle	6	2 L	00-128-020E
canister	3	5 L	00-128-050E

We protect  
your health.




Dr. Schumacher

# EXTRA enzyme-power for clean instruments

## PLURAZYME® EXTRA

Your partner for clean instruments

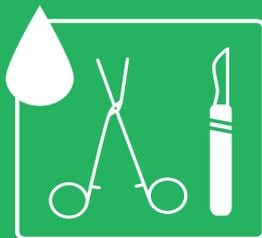

### Carefree

Although thorough cleaning is the basis for effective disinfection, the compatibility to various materials and other processing products must also meet equally high requirements. Only in this way an uncomplicated integration of a manual cleaner into the processes of modern central sterilization departments can be realized. In accordance with these requirements, we have refined PLURAZYME EXTRA.

### Gentle

The strong but equally gentle cleaning principle of PLURAZYME EXTRA is the basis for its high material compatibility. Regarding sensitive medical products, PLURAZYME EXTRA was developed to be compatible with the products for automated reprocessing THERMOSHIELD NR and THERMOSHIELD DESINFEKTANT.

### Powerful

A powerful enzyme trio of protease, amylase and lipase efficiently breaks down any organic soiling and ensures excellent cleaning in combination with the surfactant complex. A new fragrance also maintains a pleasant working environment.

Single Unit	Unit	Content	REF
dispenser bottle	12	1 L	00-138-010E
bottle	6	2 L	00-138-020E
canister	3	5 L	00-138-050E